

PartnerMic

TIPS & TRICKS




Phonak's PartnerMic picks up the voice of a speaking partner and streams the sound directly into your child's hearing aids to improve access to speech at a distance and/or in background noise.

PAIRING

The PartnerMic must first be paired to the hearing aids before use.

1. Make sure the hearing aid(s) are turned on.
2. Start the PartnerMic and hold it within 4 inches of the hearing aid(s).
3. Pairing is complete when the indicator light turns green. You will also hear a confirmation beep in the hearing aids.
4. You only have to pair the hearing aids once. They will remain paired after being turned off and on again.
5. Once paired, the PartnerMic works automatically when turned on.

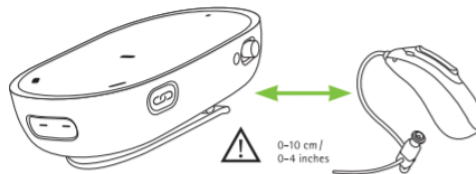

USAGE

Use the PartnerMic to send the speaker's voice directly to the hearing aids.

Consider using when listening from a distance (up to 80 feet) or in noise (in the car, at the park, during walks, or when your child is in the stroller).

1. Speaker should wear the PartnerMic approximately 8 inches below his or her mouth. The PartnerMic can be fastened to clothing or attached to the supplied lanyard.
2. Ensure the microphone is not covered by clothes or jewelry and does not bump into objects when on the lanyard.
3. The volume of the PartnerMic can be adjusted with the multi-function button on the hearing aid if activated or from the Remote Control functions within the myPhonak app.

INDICATOR LIGHTS

During operation	
Solid green 3 sec	Switching on
Blinking red	Low battery (~30min left)
Blinking blue	Pairing mode
Solid red 5 sec	Pairing failed
Staying solid red	Device error
During charging	
Blinking green	Charging
Staying solid green	Charging complete

